


TWIN / CLICK OFF

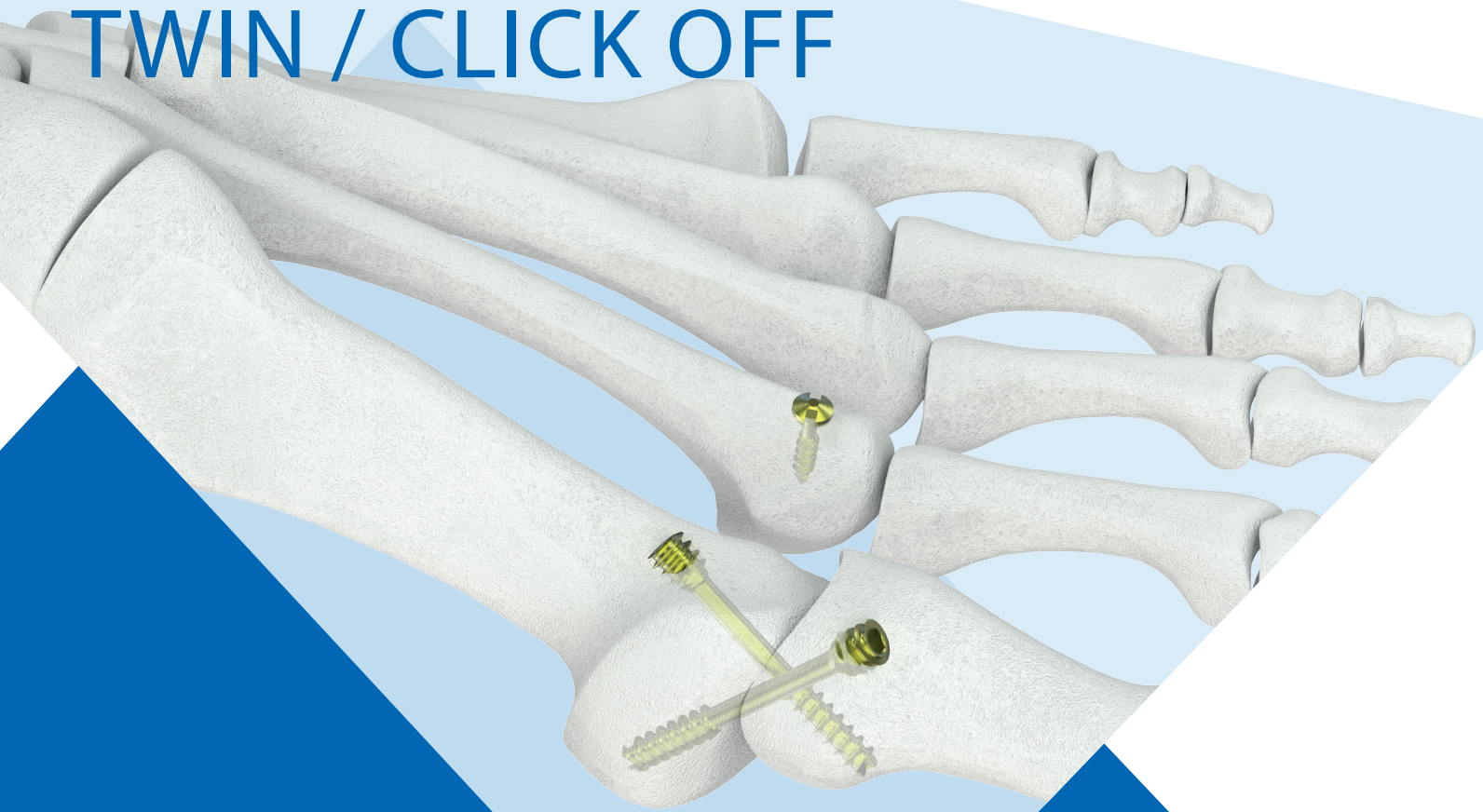

Forefoot Screw System

TWIN

Two threads, precise placement, stable compression

Cannulated compression screws in different diameters for a broad diameters for a wide range of indications. Compression is achieved by different thread pitches at the screw head and screw tip. The reverse cutting edge facilitates removal.

CLICK OFF

For fixation of small bone fragments and osteotomies in the forefoot. CLICK OFF screws are self-drilling and self-tapping. The screw is inserted directly into the bone.

TWIN Compression Screw


Compression due to
different thread pitches
at the screw head and screw tip



Reverse cutting
for easier removal when needed




Color coding GREEN / GOLD
for quick differentiation

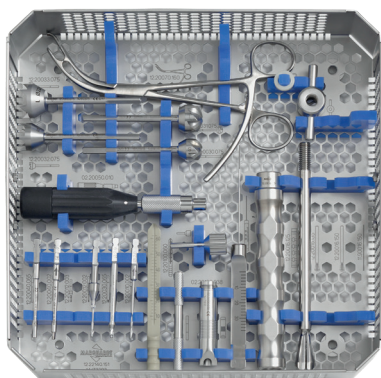

CLICK OFF Screw

self-drilling and self-cutting

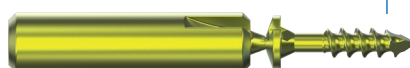

Insertion interface

TWIN / CLICK OFF Instruments



TWIN Screw Ø 2.4 / 3.2 mm, cannulated

12.03324.xxxS

Thread Length 4 - 13 mm

Length 12 - 40 mm



TWIN Screw Ø 3.0 / 4.0 mm, cannulated

12.03340.xxx(S)

Thread Length 5 - 12 mm

Length 10 - 40 mm



CLICK OFF Screw Ø 2.0

12.03402.xxx(S)

Length 10 - 17 mm



Staple Ø 1.0 mm

12.91009.xxx(S)

Angle 90°, 60°

Staple Ø 1.6 mm

12.91511.xxx(S)

Angle 90°, 60°

You can also find us on social networks:



Note:

This description is not sufficient for immediate application of the instruments. Further information is available in the detailed surgical technique. Maintenance, care and preparation of Marquardt instruments can be found in the corresponding instructions.